



The plus that gives you more

DRÄGER
FABIUS® PLUS

Every day, a new challenge



D-35161-2009

Is there a premium quality anesthesia solution out there that's economical to operate?
One that can be configured to suit our needs? One that's simple and reliable?

Take a closer look at the Dräger Fabius® plus

Hospitals around the world are constantly searching for practical, economical solutions to their everyday challenges. With more than 30,000 Fabius family members currently in use in more than 120 countries around the world, it's clear that for many hospitals the search for quality and reliability often ends with the Fabius plus. Innovative ventilation technology paired with easy operation and maintenance make the Fabius plus a natural choice for many applications.

The Fabius plus combines uncompromising flexibility and enhanced connectivity features that set it apart from so many other platforms. It gives you what you need now and still leaves room for you to grow in the future thanks to a wide range of options that let you expand your capabilities when-ever you need to. The Fabius plus is a sound investment.

PERFORMANCE FEATURES

- Supports all major ventilation modes
- High-contrast monitor with color option
- Wide-base trolley with extra stability for flexible monitor combinations
- Optional wall or ceiling versions
- Solid, spacious design with large drawers (trolley version only)
- Intelligent backup functions help maximize OR uptime
- Standardized Dräger user interface for easy and intuitive operation
- Industry standard rails for additional Dräger and third-party equipment
- CLIC adapter for use with Drägersorb Clic compatible





HIGH-PERFORMANCE VENTILATION

At its core, the Fabius plus features our powerful and accurate E-vent piston ventilator, which requires no drive gas and delivers ICU-like ventilation performance. This means the Fabius plus can handle a wide range of ventilation applications and acuity levels in pediatric and adult patients. All major ventilation modes are supported, including volume control, pressure control, pressure support (PS), and SIMV/PS. But it doesn't stop there. The Fabius plus comes complete with advanced features such as fresh gas decoupling, which maintains constant tidal volumes regardless of fresh gas flow. Dynamic compliance compensation and backup ventilation for PS modes let you ventilate with greater confidence and control. The Fabius plus can be equipped with two vaporizers for volatile anesthetic agents. The combination of its high-precision ventilator and its leak-tight breathing system design makes the Fabius plus well suited for low-flow ventilation.

In emergency situations, the quick-start feature allows manual ventilation when the device is off, on, or in standby mode for rapid support. Even in the event of a gas supply failure, automatic ventilation can be temporarily continued.

COMPREHENSIVE CONTROL

The Fabius plus features an integrated, high-contrast 6.5-inch monitor which displays all relevant ventilation parameters and a pressure curve. A color option is also available. Gas flow rates are monitored via conventional glass flow tubes.

If you've used a Dräger workstation before, then you're already familiar with Dräger's standard user interface which features our established select-adjust-confirm concept. This means you'll be up and running on the Fabius plus in just a short time.

ENHANCED ERGONOMICS AND ECONOMICS

The Fabius plus was designed to suit your work flow, not vice versa. The stable, wide-base trolley makes positioning the Fabius plus both safe and simple. If you prefer, a wall- or a ceiling-mounted version is available.

In the trolley version, three large drawers provide ample storage space for additional equipment and consumables, and a large table facilitates comfortable documentation. The design of the external compact breathing system (COSY) supports easy assembly and cleaning and can be mounted on either side to suit your individual operating room setup. Because the piston-driven E-vent ventilator requires no drive gas, no costly medical-grade gases are required for operating the ventilator.

The Fabius Plus is compatible with our CLIC adapter for use with DrägerSorb Clic which effectively prevents skin

contact with soda lime and can be exchanged even during operation. This means increasing safety for personnel and patients when replacing the sodalime and improving work processes (c.a.r.e. No. 2/2007 15)

EXPAND YOUR POSSIBILITIES: CUSTOMIZE ACCORDING TO YOUR NEEDS

Because every hospital working environment is different, the Fabius plus was designed for flexibility. Take advantage of the ample mounting space for additional monitors or IT solutions. Choose from a large selection of Dräger Infinity monitors or even from a list of Dräger-approved third-party monitors. You'll always have the option to customize the Fabius plus to suit your individual needs. If your requirements change in future, the Fabius plus can change with you. Your Dräger representative will be glad to assist you in choosing the options that are right for you.




D-35417-2009



D-35438-2009

DRÄGER FABIUS® PLUS

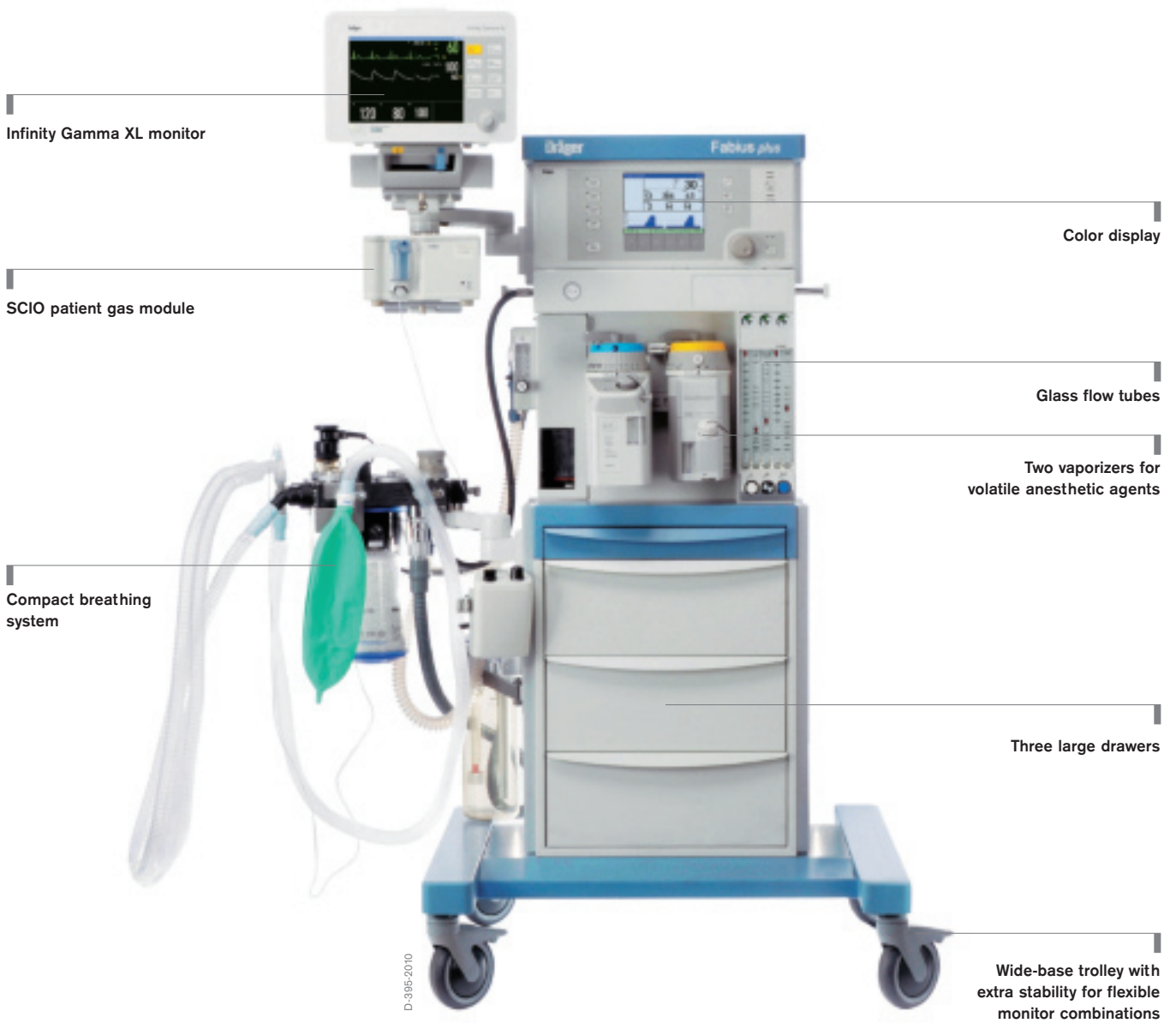

Ask your Dräger representative about other members of the Fabius family:





Now make the right choice for service: DrägerService®

DrägerService offers a wide range of flexible service options which you can select and combine to meet your individual needs – from technical support to complete, fixed-budget service contracts. Remote Service solutions offer additional potential to improve the uptime of devices.

Options and accessories

A wide array of options and accessories is available to choose from, giving you the power to create exactly the workstation you need:



D-395-2010

Color monitor option




MT-1353-2008

Breathing system with integrated heater



D-354-2010

Dräger Infinity patient monitors and mounting solutions

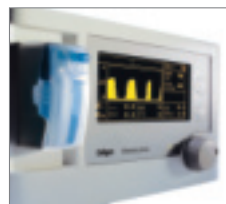

MT-24206-2009

Dräger Vamos gas monitor for volatile agents; CO₂ and SpO₂ options



D-354-2010

Auxiliary O₂ flow meter



MT-3876-2008

Flip-up side tray




MT-3890-2008

Medical-grade power strip



MT-1926-2008

Side basket

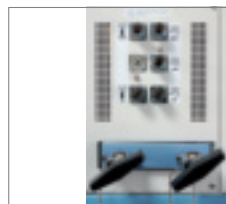

D-3-4381-2009

Pin-index mounting for reserve gas cylinders



MT-6629-2008

Dräger Remote Service Link



CORPORATE HEADQUARTERS
Drägerwerk AG & Co. KGaA
Moislinger Allee 53–55
23558 Lübeck, Germany

www.draeger.com

**REGION EUROPE CENTRAL
AND EUROPE NORTH**
Dräger Medical GmbH
Moislinger Allee 53–55
23558 Lübeck, Germany
Tel +49 451 882 0
Fax +49 451 882 2080
info@draeger.com

REGION EUROPE SOUTH
Dräger Médical S.A.S.
Parc de Haute
Technologie d'Antony 2
25, rue Georges Besse
92182 Antony Cedex, France
Tel +33 1 46 11 56 00
Fax +33 1 40 96 97 20
dlmfr-contact@draeger.com

**REGION MIDDLE EAST, AFRICA,
CENTRAL AND SOUTH AMERICA**
Dräger Medical GmbH
Branch Office Dubai
Dubai Healthcare City, P.O. Box 505108
Dubai, United Arab Emirates
Tel + 971 436 24 762
Fax + 971 436 24 761
contactuae@draeger.com

REGION ASIA / PACIFIC
Dräger Medical
South East Asia Pte Ltd
25 International Business Park
#04-27/29 German Centre
Singapore 609916, Singapore
Tel +65 6572 4388
Fax +65 6572 4399
asia.pacific@draeger.com

CANADA
Dräger Medical Canada Inc.
120 East Beaver Creek Road Suite 104
Richmond Hill Ontario L4B 4V1, Canada
Tel +1 905 763 3702
Toll-free +1 866 343 2273
Fax +1 905 763 1890
Canada.Support@draeger.com

Manufacturer:
Dräger Medical GmbH
23542 Lübeck, Germany
The quality management system at
Dräger Medical GmbH is certified
according to ISO 13485, ISO 9001
and Annex II.3 of Directive
93/42/EEC (Medical devices).